
*-Life*

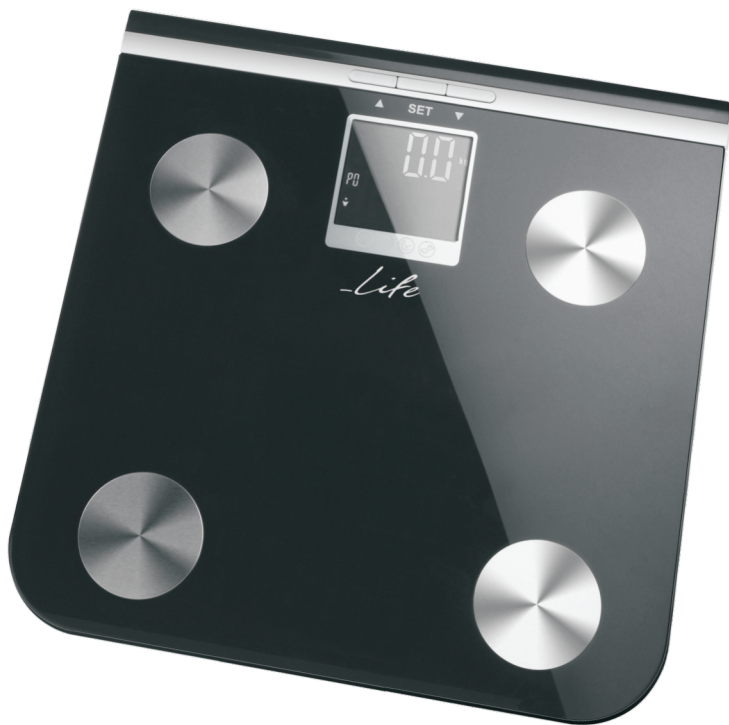**BSC-100**

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## ENGLISH

### Instruction manual

For your safety please always read the instruction manual carefully before using.

### General safety Instructions

Read the operating instructions carefully before putting the appliance into operation and keep the instructions including, the receipt and, if possible the box with the internal packing. If you give this device to other people, please also pass on the operating instruction.

- The device is designed exclusively for private and not for professional use.
- Use the device only for the intended purpose.
- Do not use a damaged device.
- Keep away from heat, direct sunlight, damp (never immerse in liquids) and sharp edges.
- Do not try to repair the appliance on your own. Always contact a qualified technician.

## Symbols in these Instructions for Use

Important information for your safety is specially marked. It is essential to comply with these instructions in order to avoid accidents and prevent damage to the machine:

### **WARNING:**

This warns you of dangers to your health and indicates possible injury risks.

### **CAUTION:**

This refers to possible hazards to the appliance or other objects.

### **NOTE:**

This highlights tips and information.

## Children

• In order to ensure your children's safety, please keep all packaging (plastic bags, boxes, polystyrene etc.) out of their reach.

### **WARNING!**

 **Caution!** Do not allow children to play with the foil as there is a danger of suffocation!

• Children should be supervised in order to ensure that they do not play with the device.

## Special Safety Instructions for this appliance

### **WARNING**

If plastic surface is moist or covered with water droplets, please wipe the surface as there is a danger of slipping.

### **CAUTION:**

- Handle the bathroom scale with care, as it is a precise measuring tool. Do not drop it or jump on it.
- Do not load the appliance with more than 150kg to avoid damaging the device.
- If the scales are carried or not used, please store them only in a horizontal position.

## Safety instructions for the batteries

### **WARNING:**

- Do not expose the battery to high temperatures or direct sunlight. Never throw batteries into the fire. Explosion danger!
- Keep batteries away from children.

Inserting the batteries

### **NOTE:**

The supplied batteries have been protected for transport with plastic foil.

1. Remove this foil before inserting the batteries into the appliance.
2. Open the battery compartment at the bottom of the

scale.

3. Insert 2 x AAA batteries. Please ensure that the polarity is correct (see the indication of the battery compartment)!
4. Close the battery compartment again.
5. If the device is not used for a prolonged period of time, please remove the batteries in order to prevent leakage of battery acid.
6. If the message "LO" appears on the display, replace the batteries.

### **CAUTION:**

- Different battery types or new and used batteries must not be used together.

## Location of controls

- 1 Contact surfaces
- 2 Display
- 3 DOWN-button ▼
- 4 UP-button ▲
- 5 SET-button

Display symbol	Meaning
<b>Kg/lb/st)</b>	Body weight
<b>P0-9</b>	User memory space
	Male
	Female
CM(FT;In)	Body height
Age	Age
FAT	Body fat display
TBW	Body water display
MUS	Muscle mass display
BONE	Bone weight display
%	Percentage-body fat/body water/muscle mass/bone weight

## Technical Specifications

Power supply: 2 x 1,5V (AAA) batteries  
Net weight of device: 2kg  
Memory capacity: 10  
Height range: 75 - 225 cm (2,5 -7,4Ft)  
Minimum weight: 2,5kg (5,5lb)  
Maximum weight: 150kg (330lb)  
Weight measurement unit: 100g increments  
Body fat measuring unit: 0,1%  
Body fat measuring range: 4,0%-50,0%  
Body water measuring range: 35,0%-85,0%  
Muscle mass measuring: 1,0%-50,0%

## Information on Body Fat measuring

### **⚠ WARNING:**

Do not use the scale if you use a pacemaker or another medical implant with electronic components. The function of the implant may be impaired by the measurement.

### **i NOTE:**

The device can only operate correctly if you stand on the contact surface with bare, dry feet!

### **Use**

The device has a very sensitive sensor. If you want to adjust the settings, press the buttons gently with your finger.

### **i NOTE:**

The scale are factory set with the units "kg" and "cm". To change the settings to "lb" or "kg" and "Ft: in", proceed as follows:

On the bottom of the scale beside the battery compartment a switch can be found. Push the switch shown below and choose between the measuring units kg/lb/st.

### **Weighting Function**

- Place the scale on a hard, flat surface (avoid carpets). An uneven floor contributes to inaccurate measurements.
- Stand on the scales with both feet. The device switches on automatically.
- Share your weight equally and wait for the scale to determine your weight.
- The digits on the display will flash before your exact weight will be shown.
- The device will automatically turn off 10 seconds after you have stepped off the scale.

### **Saving Personal Data**

#### **i NOTE:**

- The scale has default values that you can change, as described below.
- Carry out the following steps in sequence and within a few seconds. Otherwise, the device will automatically cancel input mode.
- Press the button SET. The memory space will flash on the display (USER0). Confirm the setting by pressing SET or select another memory space by pressing the UP/DOWN button. Press SET to go to the next setting.
- The symbol for «male» or «female» will flash.

Confirm the setting by pressing SET or use the UP/DOWN button to choose the symbol for your gender. Press to go the next setting.

- The height setting will flash. Select your height by pressing UP or DOWN button. Holding down one of the buttons will accelerate the process. Confirm your entry by pressing SET.
- The setting will flash. Select your age by pressing UP or DOWN button. Holding down one of the buttons will accelerate the process. Confirm your entry by pressing SET.

#### **i NOTE:**

- The device can only operate correctly if you stand on the device with bare, dry feet!
- If necessary wipe the contact surface to improve the results.
  - Due to the significant fluctuation of body fat, it is difficult to attain a standard value. On average, the optimal testing time is between 6 p.m. and 8 p.m., since the body fat remains very consistent during this period of time.
  - Press the SET button. The preset stored position date appears on the display (USER flash). Use the UP/DOWN button to select the stored position at which your date is saved.
  - Wait for the USER display and 0.0kg(0.0lb/st) to stop flashing.
- Stand on the contact surfaces (1) of the scales with both feet.
  - Firstly, your weight will flash on the display.
  - After a few seconds, the following results will be shown:
    - the body fat content in%(FAT symbol)
    - the body water content in %(TBW symbol)
    - the results analysis in the form of a bar
    - the muscle mass content in %(MUS symbol)
    - the bone mass content in %(BON symbol)
  - the present gender
  - the present height
  - the present age